ALLMOBILITY BAGNOLETTO

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BAGNOLETTO is a PVC sheet with stiffened side edges that are lifted and hooked together to form a container. It is positioned more or less in the same way as when changing the sheet with the patient on the bed. The operation is facilitated by the gripping straps.

Extend BAGNOLETTO

Position the person slowly on one side

Assembly

Place BAGNOLETTO partially under the person. Rotate the person and fully unfold BAGNOLETTO





Bathting

Use the pressure pump-tank and the hand shower. Give the pleasure of a complete bath

Dismantling

Empty the excess water. Remove the disposable sheet and disassemble BAGNOLETTO

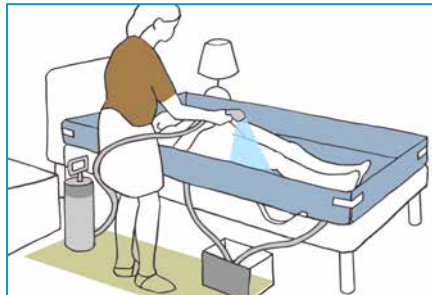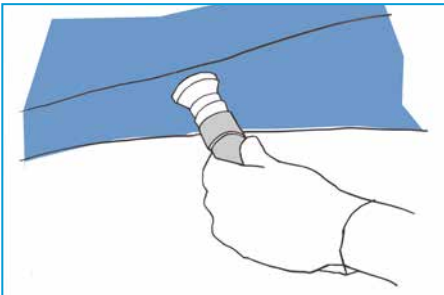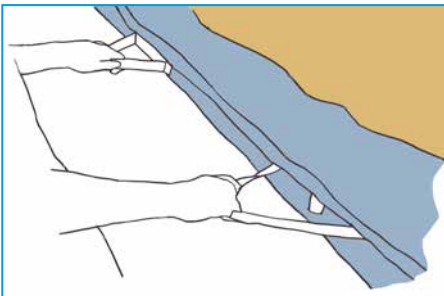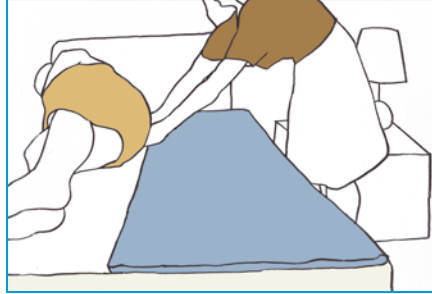The system has several water drainage points and straight or curved corrugated hoses that facilitate this operation.

The water can be brought from a tap with an extension, or the manual pressure tank (accessory) can be used. For professional use in residential or rehabilitation facilities, it is suggested to employ disposable sheets that can be used for every single wash.

Watch the video on Youtube >

ARTICLE CODE	DESCRIPTION
VQ BLBAGNOLE	BAGNOLETTO sheet 85x190x20 cm
VQ BLSERBAT	Manual pressure tank for BAGNOLETTO
VQ BLTELINI	Disposable sheets for BAGNOLETTO 10 pcs and 4 punches
VQ BLSOLOTELI	Disposable sheets for BAGNOLETTO 10 pcs

- **How to use BAGNOLETTO?**

Extend bagnoletto under your patient and use the pump-tank (which contains clean warm water) to shower the person.

- **Can it be used with any type of bed?**

Yes, BAGNOLETTO is suitable for any type of bed (even without sides).

- **What are the dimensions of the BAGNOLETTO?**

85 x 190 x 20h cm

- **Which materials are involved in the production of Bagnoletto?**

BAGNOLETTO is made of polyester fabric and class 2C PVC enamel.

- **How do I drain the waste water?**

After the bath, dirty water is drained into special containers, through the pipes provided.

- **How do I store the BAGNOLETTO?**

After each use, BAGNOLETTO must be carefully dried and sanitized with an ethanol-based solution (recommended product: Melsept spray).

- **How much water does the tank contain?**

The pump tank contains 10 litres.

- **Is it possible to add drugs to the water?**

Yes, it is possible to supply solutes and medical solvents.

- **Is it necessary to move the patient?**

No, BAGNOLETTO can be used directly on person's bed and allows you to avoid any imprudent movement.

- **How many operators are necessary to use BAGNOLETTO?**

One operator is sufficient. The sequence is short and easy to follow, as shown in the video uploaded.

- **Is it possible to use BAGNOLETTO without electricity?**

Yes, BAGNOLETTO can be used in environments without electricity.

- **Is it necessary to sanitize the aid after use?**

Disposable sheets avoid the need any further sanitizing operation, especially for healthcare facilities.